

Bovie®

Bovie Medical Corporation's electrosurgery business is now a part of



DERM 942

High Frequency Desiccator



The DERM 942 and its patent pending handpiece are designed and engineered to be the most reliable and durable desiccator available today. The unit's life is prolonged by incorporating the latest in chip technology and reducing the number of internal parts, thereby generating less internal heat.

40 watt high frequency desiccator

Adjustable in 1/10th watt increments in the first 10 watts, in monopolar and bipolar

More accurate power with Universal Power Control 100-240 VAC

Less to do and less to buy makes this the most cost effective and user friendly unit on the market

The Bovie DERM 942 is a ground reference unit which means no grounding pad is required.



BOVIE® 942 HIGH FREQUENCY DESICCATOR




STANDARD EQUIPMENT

- A942 (1) electrosurgical generator,
- A902 (1) 3-Button Autoclavable Handpiece
- A804 (5) non-sterile sharp dermal tips
- A805 (2) sterile sharp dermal tips
- A806 (5) non-sterile blunt dermal tips
- A807 (2) sterile blunt dermal tips
- A806DE (2) non-sterile DERM-Elite™ Premium blunt dermal tips
- A807DE (2) sterile DERM-Elite Premium blunt dermal tips
- A910ST (2) sterile disposable handpiece sheath
- A910 (2) non-sterile disposable handpiece sheath
- A837 (1) wall mount kit
- 110 VAC Hospital-grade power cord*
- User's guide (MC-55-239-001)*

* Not Pictured



A902

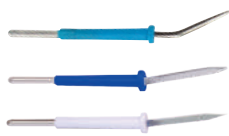

A804, A805, A806, A807,
A806DE, A807DE



A837



A910ST



A910

COMPLEMENTARY PRODUCTS †



A813
Table Top Stand †



A803
Footswitch †



A827V
Forcep Cable †



A824
Reusable Bipolar Forceps †



ES20
Ball Electrode †



ES02
Needle Electrode †



A808-T
Top Tray for A812 Stand †

A812
Mobile Stand †



† Sold Separately

PRODUCT SPECIFICATIONS

DIMENSIONS

Height: 188 mm (7.40 in.)
Width: 228 mm (8.98 in.)
Depth: 105 mm (4.13 in.)
Weight: < 2.26 kg (< 5 lbs.)

OUTPUT POWER

To 40 watts (0-10 watts in
1/10th watt increments on
both monopolar and bipolar)

OUTPUT FREQUENCY

328-492kHz

LINE VOLTAGE

A942: 100-240 VAC

LINE FREQUENCY

50 – 60Hz

WARRANTY

Four Years

